



DO LESS.

CARE MORE.



Clear Direction.  From Diagnosis to Care.

medrad® Centargo
CT Injection System

DO LESS.

Centargo minimises your injector touch time by automating manual processes.

CARE MORE.

Centargo helps you focus on what matters most in every CT scan – the patient.



With Efficiency.

Centargo is designed for a simple, reliable and automated workflow, easy daily set-up and quick change time between scans.

Centargo's daily set-up process can be completed in under 2 minutes.



The simple snap-in patient line auto-primes upon insertion and is ready for the next patient in less than 20 seconds.



With Confidence.

Less manual data entry and the intuitive user interface helps you confidently complete your daily tasks.

Touch screen within the scan room so that you can be close to your patient.



The bar code reader reduces manual data entry and provides easy traceability and access to contrast and injection details.*



* When combined with Automated documentation software



Bayer has received certification to use the Australian Made logo on the MEDRAD® Centargo Injection System, offering partners and customers confidence that they are purchasing quality innovation, supporting local jobs and helping bolster Australia's self-reliance in a critical industry.

With Centargo.

For easy set-up and patient-to-patient change over:

- › Pre-assembled 24-hour Day Set
- › Fast changeover using snap-in Patient Line

Designed with safety in mind:

- › Inlet and outlet air detection in addition to automatic air removal within reservoirs
- › Recessed connection on the Day Set
- › Two check valves in the Patient Line

Together with excellent performance:

- › Piston-based fluid delivery
- › Heat maintainers for contrast reservoirs


And modern convenience:

- › Integrated prime container
- › Dual screen configuration within scan and control room
- › Training videos on Control Room Unit
- › VirtualCare™ for remote monitoring, updates and diagnostics

With Workflow Solutions.

- › Integrated barcode reader for contrast data capture
- › Automated documentation for access to contrast and injection information in PACS
- › Protocol software tools to calculate individualised injection protocols
- › Advanced Scanner Connectivity for protocol exchange and device synchronisation*

* Requires additional ISI software

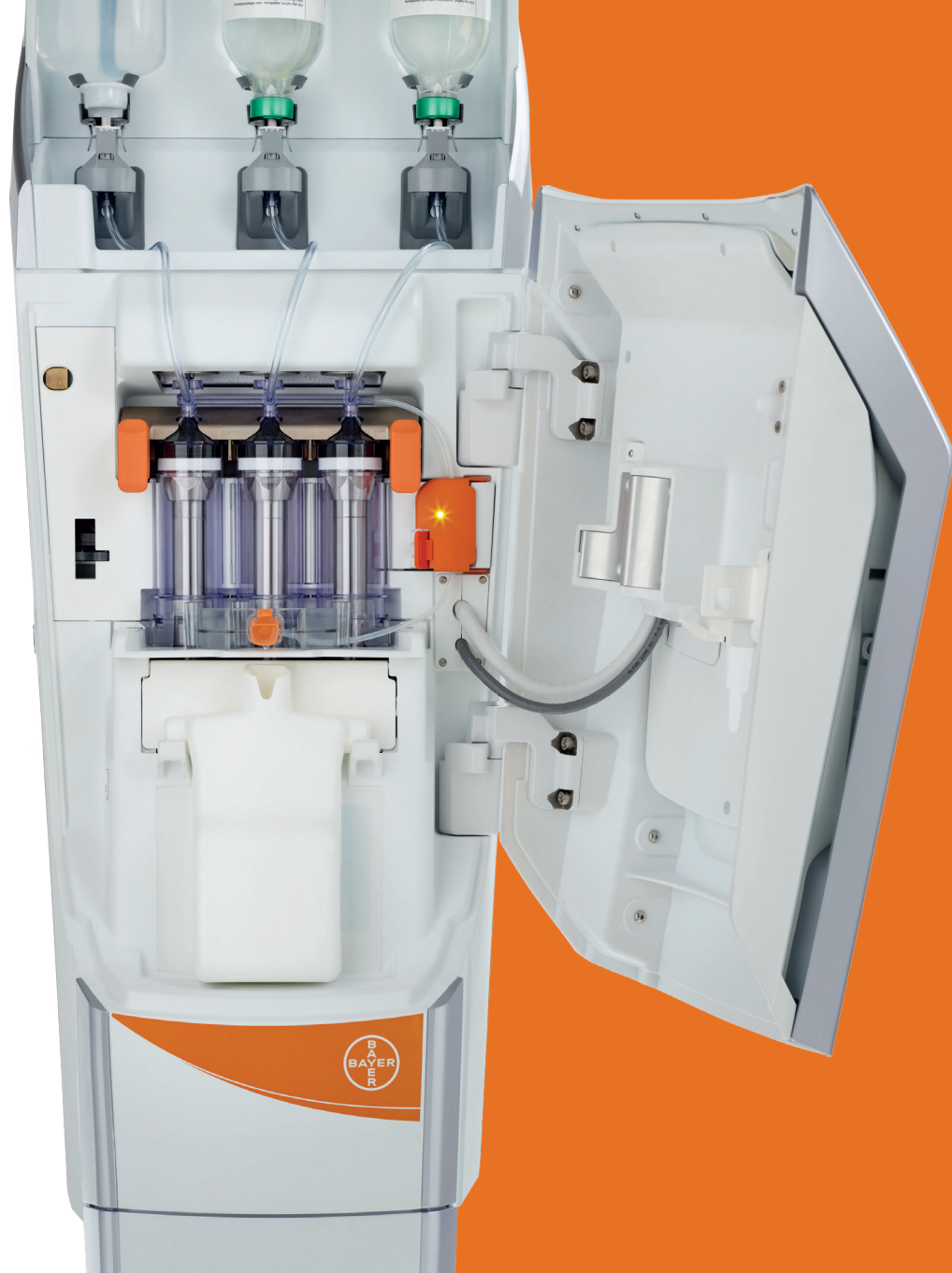

With Quality.

Centargo offers the Bayer MEDRAD® quality that you have come to trust.

- › Piston-based fluid delivery and personalised contrast media protocol software.*
- › Three inlet air sensors and one outlet air sensor as well as reservoir air detection and automatic air removal.
- › Cordless mobility through a battery powered, Wi-Fi enabled configuration.

*P3T® and Protocol Assistant Software

**Think ahead. Plan ahead. Stay ahead.
Bayer is your partner for peace
of mind and patient care.**

With Bayer.

Bayer relies on our decades of experience, comprehensive service networks, expert clinical support and dedicated sales teams to make sure that you receive the solutions you need to improve your efficiency, confidence and quality.



System Specifications

Description	Specification
Volume Range	1 mL to 200 mL for contrast and flush phases; 1 mL to 400 mL for DualFlow phases (in 1 mL increments)
Flow Rate Range	0.1 to 10 mL/sec in 0.1 mL/sec increments
DualFlow	Simultaneous injection of contrast and saline in increments of 5%
Patency Check / Test Inject	Real time flow rate adjustment and multi-phase test injections
Programmable Pressure Limit	From 50 psi to 300 psi in 1 psi increments (also displayable in kPa and kg/cm ²)
Programmable Delay	1 second to 900 seconds in 1 second increments
Pause	Maximum pause time of 20 minutes
Fluid Loading	Contrast: 50 mL – 500 mL Saline: 50 mL to 1000 mL
Fluid Reservoirs (volume)	Three 200 mL reservoirs (two contrast and one saline)
Injection Capabilities	Maximum of 6 phases per injection. Maximum of 10 injections per exam
Fluid Delivery Method	Piston-based
Day Set	Pre-assembled, labeled for 24-hour use, designed for set-up in less than 2 minutes
Patient Line	250 cm, single hand design, two check valves, auto-primers upon insertion
Displays	Two independent, Wi-Fi enabled, full colour, touch screen units with modern graphical user interfaces located in the scan room and control room
Air Detection	Three inlet air sensors and one outlet air sensor includes reservoir air detection and automatic air removal
Battery Performance	Up to 16 hours or up to 60 injections
Integrated Prime Container (vol)	680 mL
Injector Dimensions (WxDxH)	76.2 x 68.6 x 167.6 cm
Injector Weight (approx.)	76 kg (Pedestal configuration)
Available Scanner Interface	ISI2 and Connect.CT
Smart Protocols	P3T® algorithms calculate individualised injection protocols for cardiac CT, CT pulmonary angiography and CT imaging of abdominal organs Protocol Assistant provides basic protocol calculations for a range of procedure types
eGFR Calculator	Available on Control Room and Scan Room Units

Ordering Information

Injector Configurations

- Centargo Pedestal System with Batteries **CENT-SYS-BAT**
- Centargo OCS System **CENT-SYS-OCS**

Disposables

- Centargo Day Set (24-hour use) **CENT-DS**
- Centargo Patient Line **CENT-PL**
- Centargo Replacement Spikes (optional) **CENT-RS**

Smart Protocols

- P3T Personalised Protocols **CENT-P3T**
- PAT Protocol Assistant **CENT-PAT**

Workflow Solutions

- Automated Documentation **CENT-AUTO-DOC**
- Smart Protocols **CENT-SMART**
- Insights **CENT-INSIGHTS**

Processing Productivity

- RIS Interface **CENT-RIS**
- Dictation Interface **CENT-SR**
- Radimetrics Interface **CENT-REP**

Scanner Connectivity

- ISI2 **CENT-ISI2**
- Connect.CT **CENT-ISI-CCT**



More information on
radiology.bayer.com.au



Manufacturer

Imaxeon Pty Ltd
Unit 1, 38 – 46 South Street
Rydalmere, NSW 2116
Australia
Phone: 02 8845 4999
Fax: 02 8845 4998
Imaxeon is a Bayer subsidiary

Australian Sponsor

Imaxeon Pty Ltd
Unit 1, 38 – 46 South Street
Rydalmere, NSW 2116
Australia
Phone: 02 8845 4999
Fax: 02 8845 4936
Customer Service: 1800 633 723



Bayer reserves the right to modify the specifications and features described herein or to discontinue any product or service identified in this publication at any time without prior notice or obligation. Please contact your authorised Bayer representative for the most current information.

All patient data that appear in this document are fictitious. No actual patient information is shown.

Bayer, the Bayer Cross, MEDRAD, MEDRAD Centargo, Centargo, P3T, VirtualCare and Radimetrics are trademarks owned by and/or registered to Bayer in the U.S. and/or other countries. Other trademarks and company names mentioned herein are properties of their respective owners and are used herein solely for informational purposes. No relationship or endorsement should be inferred or implied.

© 2023 Bayer. This material may not be reproduced, displayed, modified or distributed without the express prior written consent of Bayer.